

Important: This product is to be used ONLY under the supervision of a medical professional.



Use the PTQ series PegAssist™ with:

# PegAssist™ PTQ Series Off-loading Insole

## PegAssist™ Off-loading Insoles

The DARCO PegAssist™ System features a multi-indication removable-peg chassis that effectively off-loads the plantar aspect of the foot after surgery or when wounds are present.

- Removable pegs allow for localized off-loading of wounds and ulcerations of the foot.
- Pressure is reduced by as much as 60% while allowing patient to remain ambulatory.
- Each insole has a 15mm combination Plastazote® and multi-foam base with a 3mm PORON® Cushioning cover.
- Included stabilizer board prevents adjacent pegs from shifting.

### Contents:

- PegAssist™ Insole
- Stabilizer board
- Red marker



### Video demonstration:


Proper assembly and application of the PegAssist™ Off-loading Insole

### Scan QR code or visit:

[darcointernational.com/pegassist](http://darcointernational.com/pegassist)

### PQ Series PegAssist™ Sizing:

Size:	Product Code:			
	Women:	Men:		
S	4 - 6	6 - 8	PTQW1	PTQM1
M	6.5 - 8	8.5 - 10	PTQW2	PTQM2
L	8.5 - 10	10.5 - 12	PTQW3	PTQM3
XL	-	12.5 - 14	-	PTQM4

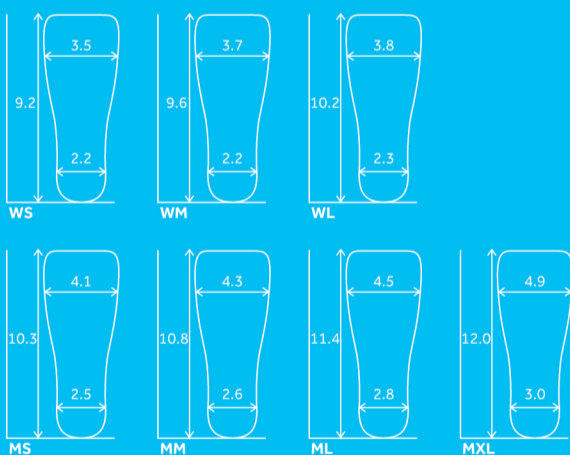

Important: This product is to be used ONLY under the supervision of a medical professional.



**Manufactured for:**  
DARCO International, Inc.  
810 Memorial Blvd.  
Huntington, WV 25701, USA  
Toll-free: 800.999.8866  
Phone: 304.522.4883  
[www.darcointernational.com](http://www.darcointernational.com)



DARCO (Europe) GmbH  
Gewerbegebiet 18  
82399 Raisting, Germany  
Tel. +49 (0) 88 07. 92 28-0  
[www.darco.de](http://www.darco.de)

Made in China

